

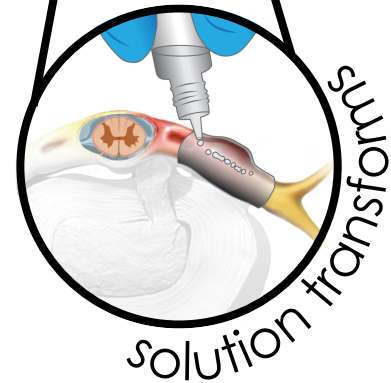


Designed for the management of **nerve root tethering** in decompression procedures

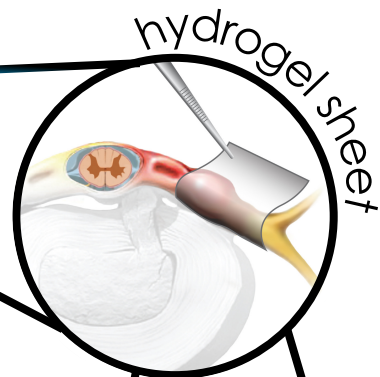

- Provides a hydrogel **interface** between the nerve root and surrounding tissues
- Advanced solution to **unwanted soft tissue attachment**
- **Not a biologic** and contains no tissue
- **Made of hyaluronic acid and alginate**
- **Tissue-adherent** across entire surface, reducing likelihood of device migration

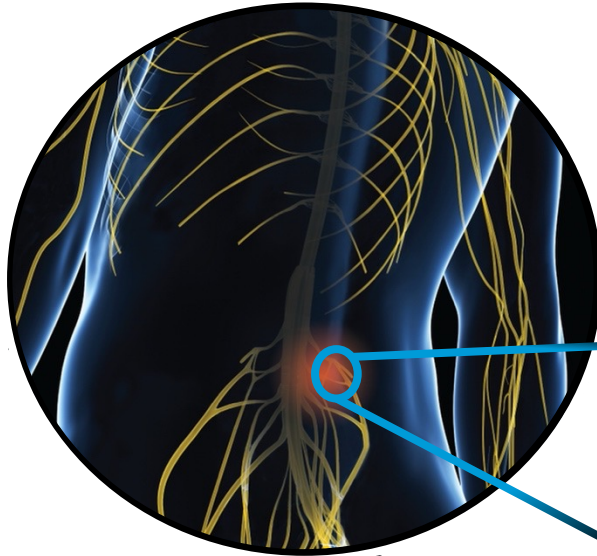


VersaWrap® is easy to apply

<p>ULTRATHIN</p> <p>The VersaWrap sheet is an ultrathin plant-based hydrogel consisting of biopolymers hyaluronic acid and alginate</p>		 <p>VersaWrap sheet placed on nerve root</p>
<p>SUTURELESS</p> <p>The VersaWrap solution is added to the sheet, rendering a tissue-adherent hydrogel layer</p>		 <p>VersaWrap solution added</p>

VersaWrap® is a bioresorbable hydrogel

VersaWrap® bioresorbs via hydrolysis and metabolism and is not remodeled or replaced, allowing tissues to return to a pre-injured state^{1,2}

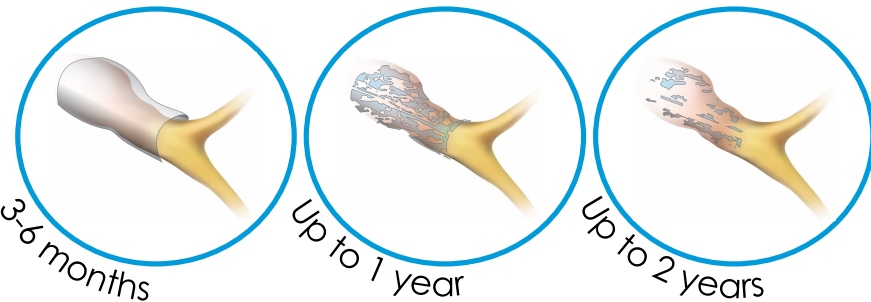

VersaWrap® is a medical device

VersaWrap® is indicated for the management of peripheral nerve injuries in which there has been no substantial loss of nerve tissue. VersaWrap® is also indicated for the management and protection of tendon injuries in which there has been no substantial loss of tendon tissue. The device may also be used in the management and protection of surrounding tissues such as skeletal muscle and ligament

VersaWrap® Ordering Information

Catalog Number: VTP-2201, VTP-1201

For order inquiries: 800.206.5586 or sales@alafairbiosciences.com

Always refer to the package insert, product label, and/or user instructions before using any Alafair product. Surgeons must rely on their own clinical judgement when deciding which treatments and procedures to use with patients. Products may not be available in all markets. Product availability is subject to the regulatory or medical practices that govern individual markets. Please contact your Alafair representative if you have questions about the availability of Alafair products in your area.

VersaWrap® is a registered trademark of Alafair Biosciences, Inc.

Copyright 2023 Alafair Biosciences, Inc.
Printed in the USA

1 TR-011, Biocompatibility, Data on file at Alafair
2 TR-018, Data on file at Alafair

MKT-04052 Rev 03